

## STAIR-TRAC

portable wheelchair lift



The Garaventa Stair-Trac is an easy to use, portable wheelchair lift designed to attach under most manual wheelchairs. Stair-Trac's powerful motor allows a small attendant to safely transport a larger passenger up and down stairways. Stair-Trac is suitable for indoor or outdoor use in public buildings and private homes. An affordable and reliable portable access solution.



## Garaventa Super-Trac

- Portable inclined platform lift
- Suitable for indoor and outdoor use
- Accommodates all types of wheelchairs

## Specifications

Load Capacity.....130kg/285lbs  
 Power Capacity.....650 stairs  
 Maximum Speed:  
     Up Stairs.....6.5m (21.2ft)/min  
     Down Stairs.....8.8m (28.9ft)/min  
 Motor.....12 VDC  
 Batteries.....2 x 6V - 20AH  
 Max. Charge Rate.....4 amp  
 Optimum Stair Angle.....35°

## Dimensions

Height.....922mm/36.3in  
 Length.....1468mm/57.8in  
 Width.....636mm/25.0in

## Weight

Drive Frame.....39.5kg/87lbs  
 Support Handle.....14kg/31lbs  
 Total Weight.....53.5kg/118lbs

## Garaventa Lift

PO Box 1769,  
 Blaine, WA 98231-1769  
 7505 134A Street, Surrey, BC  
 Canada V3W 7B3

Phone: +1 604 594-0422  
 Toll Free: 1-800-663-6556  
 Fax: +1 604 594-9915

Email: [productinfo@garaventa.ca](mailto:productinfo@garaventa.ca)  
 Web site: [www.garaventa.ca](http://www.garaventa.ca)  
 Global Portal: [www.garaventlift.com](http://www.garaventlift.com)

©2008 Garaventa Lift. As we are continuously improving our products, specifications outlined in this brochure are subject to change without notice.

Printed in Canada  
 15803-F-PB

# STAIR-TRAC - Portable Wheelchair Lift

Need Design Assistance?

Call your local Garaventa Lift representative or our Design Hotline at 1-800-663-6556

## Features


- Attaches to most standard wheelchairs with push handle (adult/child)
- Powered to carry over 250 lbs. up approximately 30 floors
- Use on virtually any indoor or outdoor stairway
- Stable and self-supporting
- Durable tracks grip stairs securely
- Detachable support handle for compact storage and easy transport
- Quick-release safety belt
- Key switch restricts use to authorized persons
- Rechargeable batteries and charger
- Portable inclined platform lift
- Suitable for indoor and outdoor use
- Accommodates all types of wheelchairs

## Benefits

- Economical compared to stairway lifts and elevators
- It is portable so it can be used on various stairways within a building
- Can be used when electrical power fails or is not available
- Easily stores in a secure area when not in use
- Requires minimal maintenance
- Available immediately
- Can be used for building evacuation

## Clearances

Wheelchair loading/unloading requires a minimum straight line clearance of approximately 3m/10ft. The diagram below shows typical clearances required on turnback stairways. Dimensions vary with stairway width and operator experience. Consult Garaventa or your local representative for clearance verification.

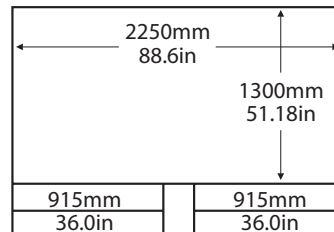

## What our customers say

*"The Garaventa Stair-Trac has allowed our District to meet its needs at a fraction of the cost of installing an elevator in our building."*  
 - Batavia City Schools, Batavia, NY

*"It has proven to be dependable, safe, and functional for our school district."*  
 - Hillcrest Elementary, American Falls, ID



Attach the Stair-Trac to the wheelchair, recline the passenger and move them to the stairway.



Ascend or descend stairs with the press of a button. The passenger sits in a comfortable, upright position.



Stair-Trac will climb virtually any stair tread material.



[www.garaventa.ca](http://www.garaventa.ca)

Authorized Garaventa Lift Representative

