

STAIR-TRAC

portable wheelchair lift



The Garaventa Stair-Trac is an easy to use, portable wheelchair lift designed to attach under most manual wheelchairs. Stair-Trac's powerful motor allows a small attendant to safely transport a larger passenger up and down stairways. Stair-Trac is suitable for indoor or outdoor use in public buildings and private homes. An affordable and reliable portable access solution.



STAIR-TRAC - Portable Wheelchair Lift

Need Assistance?

Call your local Garaventa Lift representative or our Hotline at 1-800-663-6556

Features

- Attaches to most standard wheelchairs with push handle (adult/child)
- Powered to carry over 250 lbs. up approximately 30 floors
- Use on virtually any indoor or outdoor stairway
- Stable and self-supporting
- Durable tracks grip stairs securely
- Detachable support handle for compact storage and easy transport
- Quick-release safety belt
- Key switch restricts use to authorized persons
- Rechargeable batteries and charger
- Portable inclined platform lift
- Suitable for indoor and outdoor use
- Accommodates all types of wheelchairs



Attach the Stair-Trac to the wheelchair, recline the passenger and move them to the stairway.

Garaventa Super-Trac

- Portable inclined platform lift
- Suitable for indoor and outdoor use
- Accommodates all types of wheelchairs

Specifications

Load Capacity.....130kg/285lbs
 Power Capacity650 stairs
 Maximum Speed:
 Up Stairs.....6.5m (21.2ft)/min
 Down Stairs8.8m (28.9ft)/min
 Motor12 VDC
 Batteries.....2 x 6V - 20AH
 Max. Charge Rate4 amp
 Optimum Stair Angle.....35°

Dimensions

Height.....922mm/36.3in
 Length1468mm/57.8in
 Width636mm/25.0in

Weight

Drive Frame.....39.5kg/87lbs
 Support Handle.....14kg/31lbs
 Total Weight53.5kg/118lbs

Benefits


- Economical compared to stairway lifts and elevators
- It is portable so it can be used on various stairways within a building
- Can be used when electrical power fails or is not available
- Easily stores in a secure area when not in use
- Requires minimal maintenance
- Available immediately
- Can be used for building evacuation



Ascend or descend stairs with the press of a button. The passenger sits in a comfortable, upright position.

Clearances

Wheelchair loading/unloading requires a minimum straight line clearance of approximately 3m/10ft. The diagram below shows typical clearances required on turnback stairways. Dimensions vary with stairway width and operator experience. Consult Garaventa or your local representative for clearance verification.

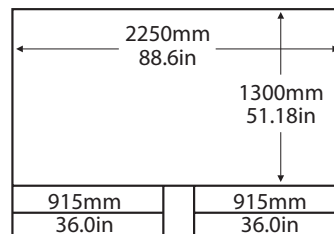

Stair-Trac will climb virtually any stair tread material.

What our customers say

"The Garaventa Stair-Trac has allowed our District to meet its needs at a fraction of the cost of installing an elevator in our building."
 - Batavia City Schools, Batavia, NY

"It has proven to be dependable, safe, and functional for our school district."
 - Hillcrest Elementary, American Falls, ID

Garaventa Lift

Toll Free: 1.800.663.6556
 Email: productinfo@garavalift.com
 Web site: www.garavalift.com

Global Portal: www.garavaliftgroup.com

© Garaventa Lift. As we are continuously improving our products, specifications outlined in this brochure are subject to change without notice.

Printed in Canada
 15803-G-PB



www.garavalift.com

Authorized Garaventa Lift Representative

